


NEW!
St. Jude Medical
CPS™ model available.

Pressure Products, Inc.

H. Lee **SafeSheath®** Sealing Adapter

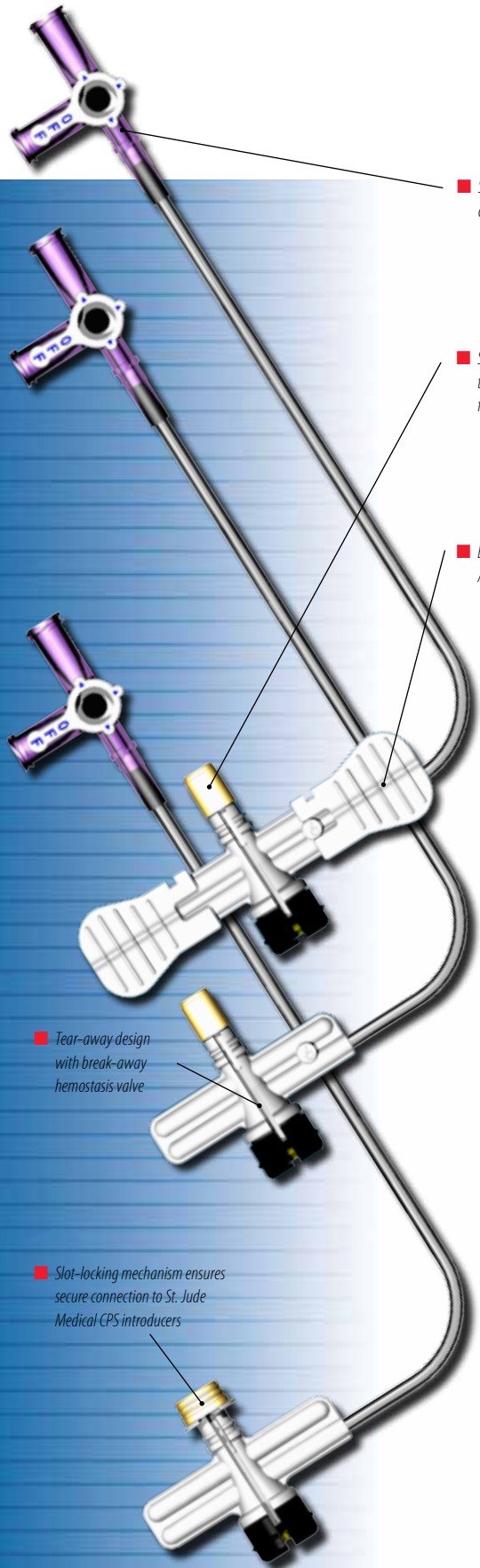

■ Side port for infusion and contrast injection

■ Securely attaches to introducer using the press-fit adapter, providing a fluid- and air-tight connection.

■ Extended Wing version for Medtronic Attain for easier hub separation

■ Tear-away design with break-away hemostasis valve


■ Slot-locking mechanism ensures secure connection to St. Jude Medical CPS introducers

Provides the convenience of a splittable, hemostatic valve when using a Medtronic Attain™ or St. Jude Medical CPS™ introducer.

The **SafeSheath Sealing Adapter** allows you to combine the advantages of the SafeSheath system while employing a Medtronic Attain™ or St. Jude Medical CPS™ introducer.

Features:

- Unique press-fit connector provides a fluid- and air-tight seal
- Innovative tear-away design with break away hemostasis valve
- Hemostasis when single or multiple wires are within it, as well as when leads need to be introduced or exchanged
- Reduced risk of blood loss and air embolism
- A side port for “keep open” infusion

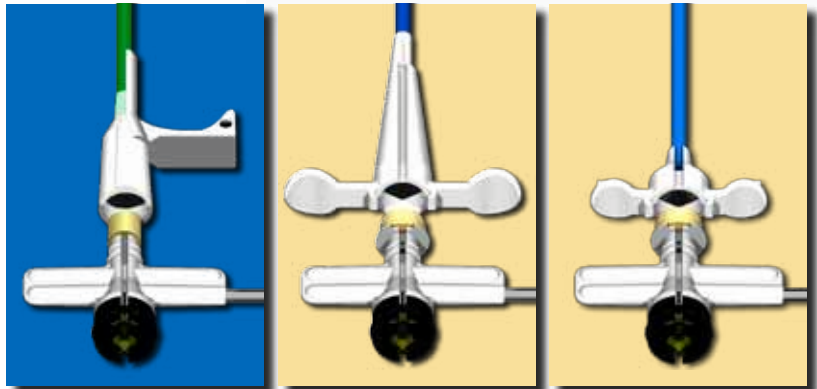

SafeSheath Sealing Adapter model for Medtronic Attain™ introducer. Also available with extended wings.

SafeSheath Sealing Adapter model for St. Jude Medical CPS™ braided and non-braided introducers.

Model	Description	Size
<i>For Medtronic Attain™ Introducer</i>		
SS-SA-09	SafeSheath Sealing Adapter	9F
SS-SA-EW-09	Extended Wing SafeSheath Sealing Adapter	9F
<i>For St. Jude Medical CPS™ Introducer</i>		
410195	SafeSheath Sealing Adapter	9F

The H. Lee SafeSheath™ hemostatic splittable sheath is protected by Patent Insurance through 2011. USA Patents 5,125,904, 5,312,355 and 5,409,463. Other patents pending. Canadian Patent 2,112,394. European Patent EPOS93685. SafeSheath is a registered trademark of Pressure Products, Inc. Attain is a trademark of Medtronic, Inc. CPS is trademark of St. Jude Medical.

P
R
O
D
U
C
T
I
N
F
O
R
M
A
T
I
O
N