


H. Lee SafeSheath® D-Pro

DIALYSIS INTRODUCER SERIES

Hemostatic, peel-away, dialysis introducer for venous access. Standard and side-port models available in 14, 15 and 16 F sizes.

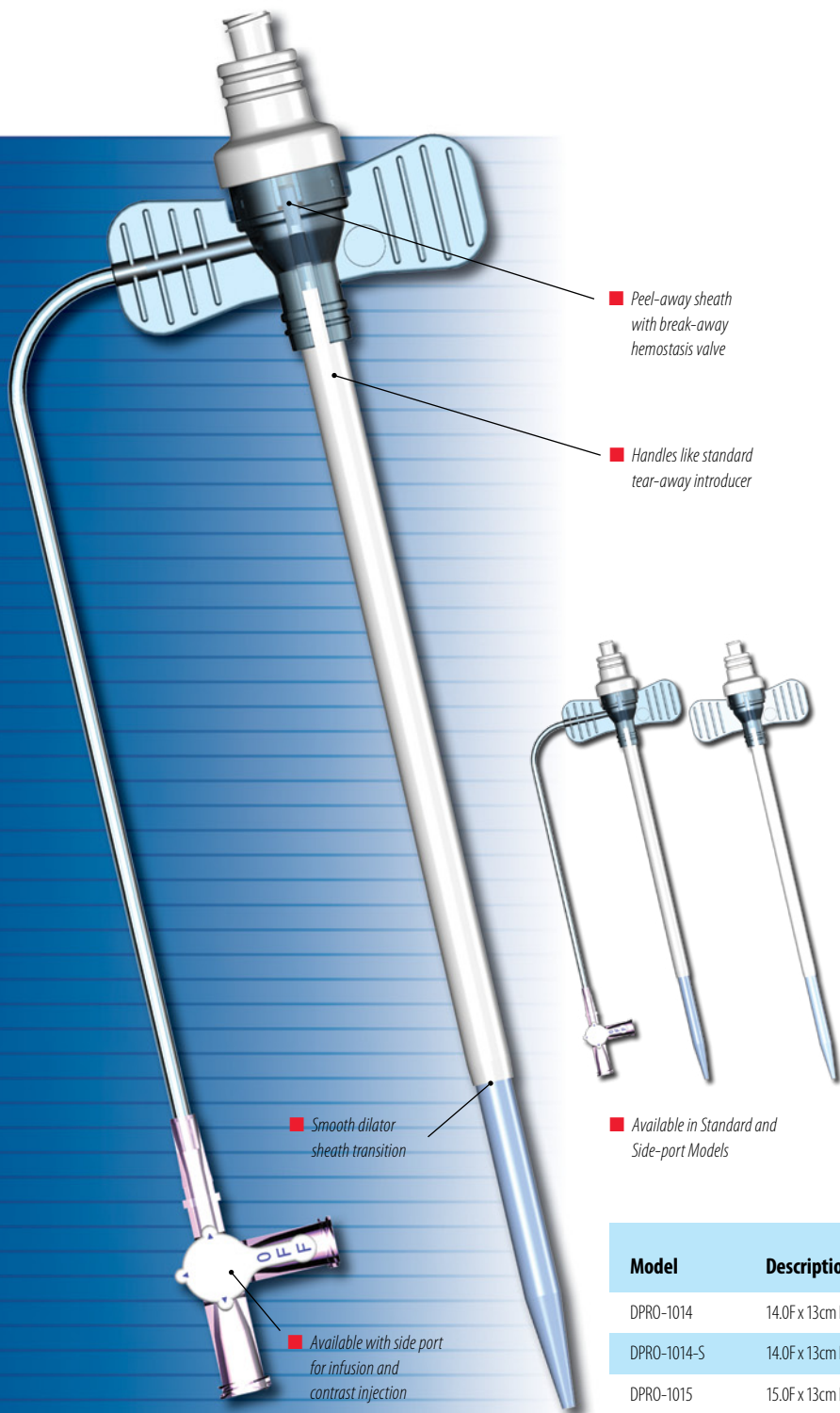

Introducing Confidence

Building on a decade of innovation as the acknowledged leader in cardiovascular hemostatic-valved peel-away sheath technology, Pressure Products proudly announces the introduction of the **SafeSheath® D-Pro** hemostatic peel-away dialysis introducer system.

SafeSheath D-Pro provides the safety and security of hemostatic access for dialysis catheter introduction by combining the convenience of a tear-away introducer with a hemostatic valve, preventing back-bleeding and greatly reducing the possibility of air embolism.

Available in standard and side port models, SafeSheath D-Pro is the intelligent choice for your dialysis access needs.

Features:

- Innovative peel-away sheath with break-away hemostasis valve
- Reduces risk of blood loss and air embolism during venous access
- Smooth dilator/sheath transition to reduce vein entrance forces
- Available with optional side port for “keep open” infusion and contrast injection
- Handles like standard tear-away introducer
- Available in 14, 15, and 16 French; additional custom sizes available

Model	Description	Dilator Length
DPRO-1014	14.0F x 13cm Peel-away hemostatic-valved dialysis sheath	20cm
DPRO-1014-S	14.0F x 13cm Peel-away hemostatic-valved dialysis sheath w/ side-port	20cm
DPRO-1015	15.0F x 13cm Peel-away hemostatic-valved dialysis sheath	20cm
DPRO-1015-S	15.0F x 13cm Peel-away hemostatic-valved dialysis sheath w/ side-port	20cm
DPRO-1016	16.0F x 13cm Peel-away hemostatic-valved dialysis sheath	20cm
DPRO-1016-S	16.0F x 13cm Peel-away hemostatic-valved dialysis sheath w/ side-port	20cm

Each set comes with a sheath and dilator packaged in a pouch, 5 to a box. The H. Lee SafeSheath® D-Pro hemostatic splittable sheath is protected by Patent Insurance through 2011. USA Patents 5,125,904, 5,312,355 and 5,409,463. Canadian Patent 2,112,394. European Patent EP0593685. SafeSheath is a registered trademark of Pressure Products, Inc.