


# SafeSheath® Ultra Lite Kit

Hemostatic Peel-away Introducer System for Vascular Access

The specially engineered membrane of the Ultra Lite's hemostatic valve provides minimal resistance during lead introduction and positioning, while reducing back bleeding and the possibility of air embolism.

*with low resistance valve, ergonomically-designed, easy-splitting hub, extruded score line sheath, and positive-locking dilator*

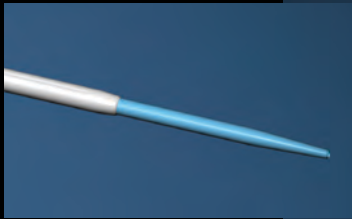

The H. Lee SafeSheath® Ultra family of hemostatic introducers continues to revolutionize venous access for device lead placement, setting new standards for performance and features while facilitating better outcomes.

#### FEATURES

- ▶ Specially designed membrane allows effortless insertion, positioning, and removal of leads and catheters
- ▶ Extruded score line peel-away sheath for clean, even peeling during removal
- ▶ Ergonomic, easy-splitting break-away hub
- ▶ Hemostatic valve permits multiple lead insertions without further interventions
- ▶ Standard 13cm length with a wide range of French sizes available
- ▶ Positive locking dilator connector
- ▶ Reduced risk of blood loss and air embolism
- ▶ French size printed on handle and indicated by color-coded introducer cap

#### PRODUCT DESCRIPTION

The SafeSheath Ultra family of introducers combines a tear-away sheath and hemostatic valve to reduce back bleeding and the possibility of air embolism.

The SafeSheath Ultra Lite is introduced into the subclavian vein in the standard manner. The hemostatic valve of the 13cm Ultra Lite seals the device during procedures.

Incorporating a specially engineered valve membrane, the Ultra Lite provides minimal insertion forces during lead introduction and positioning.



The easy splitting, ergonomically shaped break-away hub greatly simplifies sheath removal after lead or catheter placement. The extruded score line construction ensures a clean, even peel every time. A new locking dilator connector assures a firm attachment to the sheath.

Available in a wide-range of French sizes,

SafeSheath Ultra Lite is the intelligent choice for your vascular access needs.


# SafeSheath® Ultra Lite Kit

## Hemostatic Peel-away Introducer System for Vascular Access

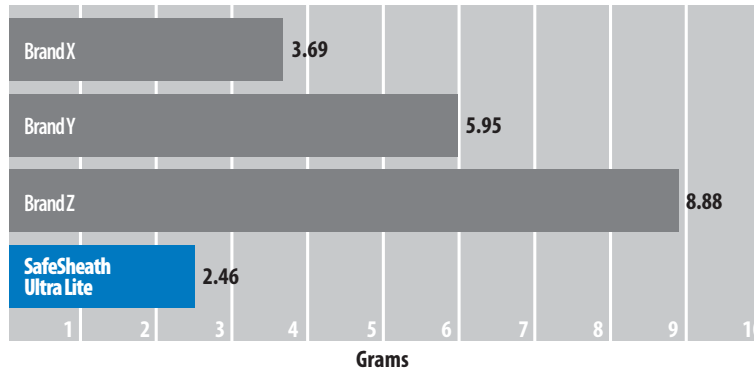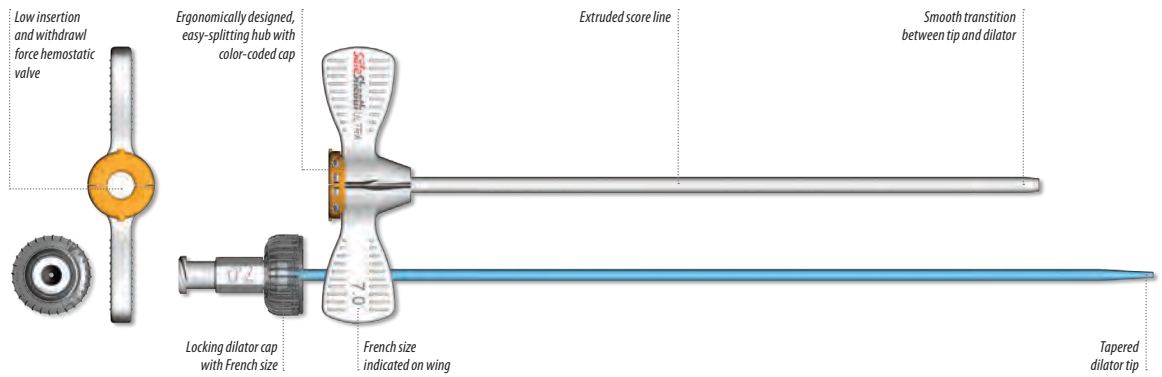

#### INSERTION FORCE REQUIRED FOR 6.6F SILICONE CATHETER

The specially developed membrane of the SafeSheath Ultra Lite valve generates extremely low insertion forces.

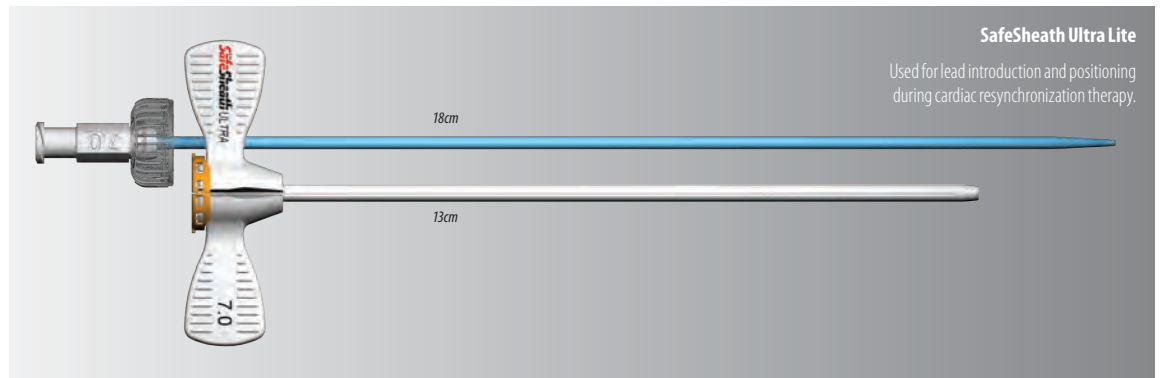

#### SafeSheath Ultra Lite Tear-away valved hemostatic sheath system

Model	Size	Hub Color	Length	Sheath I.D.	Tip I.D.	Dilator Length
SUPK5	5.0 French	Gray	13cm	.078"	.073"	18cm
SUPK6	6.0 French	Green	13cm	.090"	.085"	18cm
SUPK7	7.0 French	Orange	13cm	.102"	.095"	18cm
SUPK8	8.0 French	Blue	13cm	.113"	.108"	18cm
SUPK85	8.5 French	Light Blue	13cm	.120"	.115"	18cm
SUPK9	9.0 French	Black	13cm	.126"	.121"	18cm
SUPK10	10.0 French	Fuschia	13cm	.139"	.134"	18cm
SUPK11	11.0 French	Yellow	13cm	.152"	.147"	18cm

Guide Wire for SafeSheath Ultra Lite is 50cm long with O.D. of .038"

#### ORDERING INFORMATION

Each SafeSheath Ultra Lite Introducer System kit includes:

- ▶ One tear-away sheath
- ▶ One dilator
- ▶ One 18g needle
- ▶ One 10cc syringe
- ▶ One 50cm guidewire

Five (5) kits per box

The H. Lee SafeSheath®, hemostatic splittable sheath is protected by Patent Insurance. European Patent EP0593685.

SafeSheath is a registered trademark of Pressure Products, Inc.



#### Pressure Products, Inc.

Customer Service  
1861 N. Gaffey Street Suite B San Pedro, CA 90731  
(800) 600-4973 (310) 547-4973 Fax (310) 547-4760

www.pressure-products.com