


# SafeSept<sup>®</sup>

## Transseptal Guidewire

*Specially designed for improved safety to cross the interatrial septum at the fossa ovalis when supported by a transseptal introducer.*

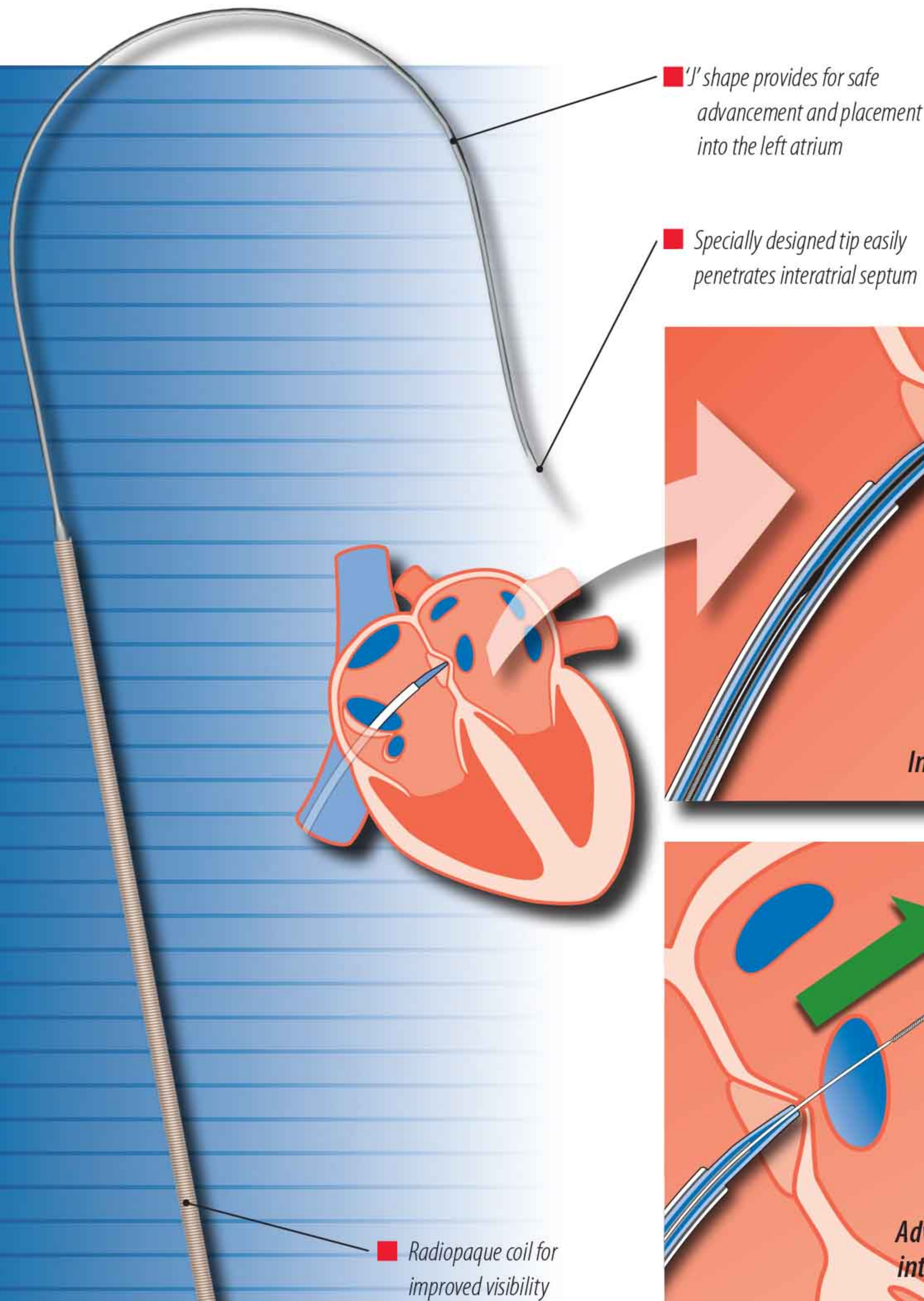

### Cross with Confidence

Introducing the 0.018" by 140cm SafeSept<sup>®</sup> Transseptal Guidewire, the next significant advance in left atrial access since the original SafeSept<sup>®</sup> Transseptal Guidewire.



SafeSept is a 140cm long, 0.018" diameter nitinol guidewire specifically designed for transseptal puncture. After the transseptal dilator has "tenting" the fossa ovalis, effortless advancement of the SafeSept tip perforates the membranous fossa. Unsupported by the needle and dilator, the tip of the wire assumes a 'J' shape, rendering it incapable of further tissue penetration.

A radiopaque coil along the shaft allows for fluoroscopic visualization of the wire within the left atrium while proximal marker bands help determine the approximate location of the SafeSept tip in relation to the tip of the needle. The transseptal needle, dilator and sheath are then advanced "over the wire" into the left atrium with virtually no possibility of aortic or pericardial perforation.

The SafeSept Transseptal Guidewire is part of a complete line of implantable device therapy products and accessories from Pressure Products.

#### Features:

- Radiopaque coil makes the small diameter 0.018 inch wire easily visible fluoroscopically
- Specially designed tip configuration enables effortless penetration of the fossa
- 'J' shape allows safe advancement to assure proper wire location within left atrium
- Transseptal needle and dilator are rendered atraumatic when tracking over the SafeSept
- Proximal marker bands help determine SafeSept tip location relative to transseptal needle tip during the procedure
- Compatible with all currently available adult transseptal introducer systems

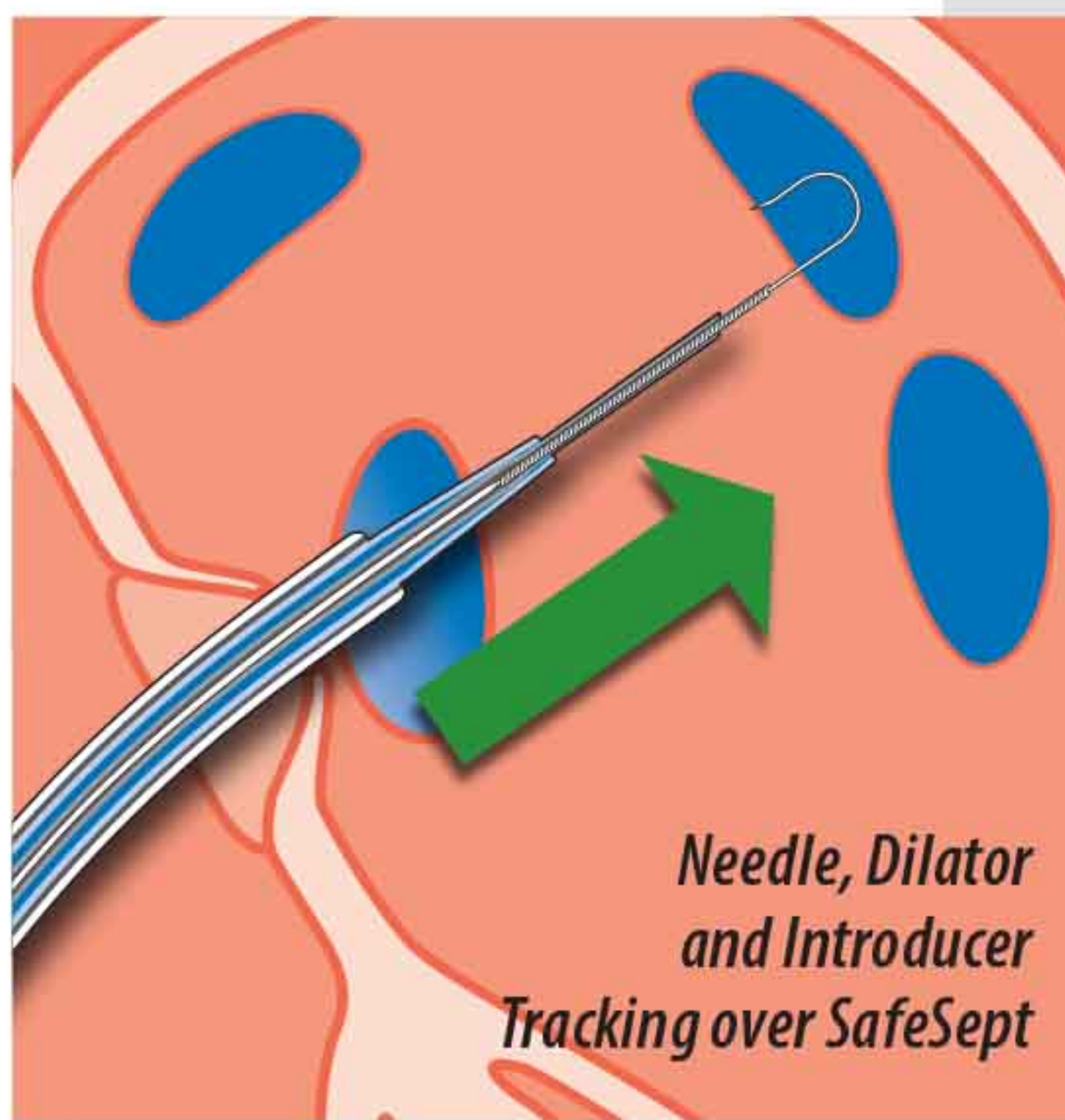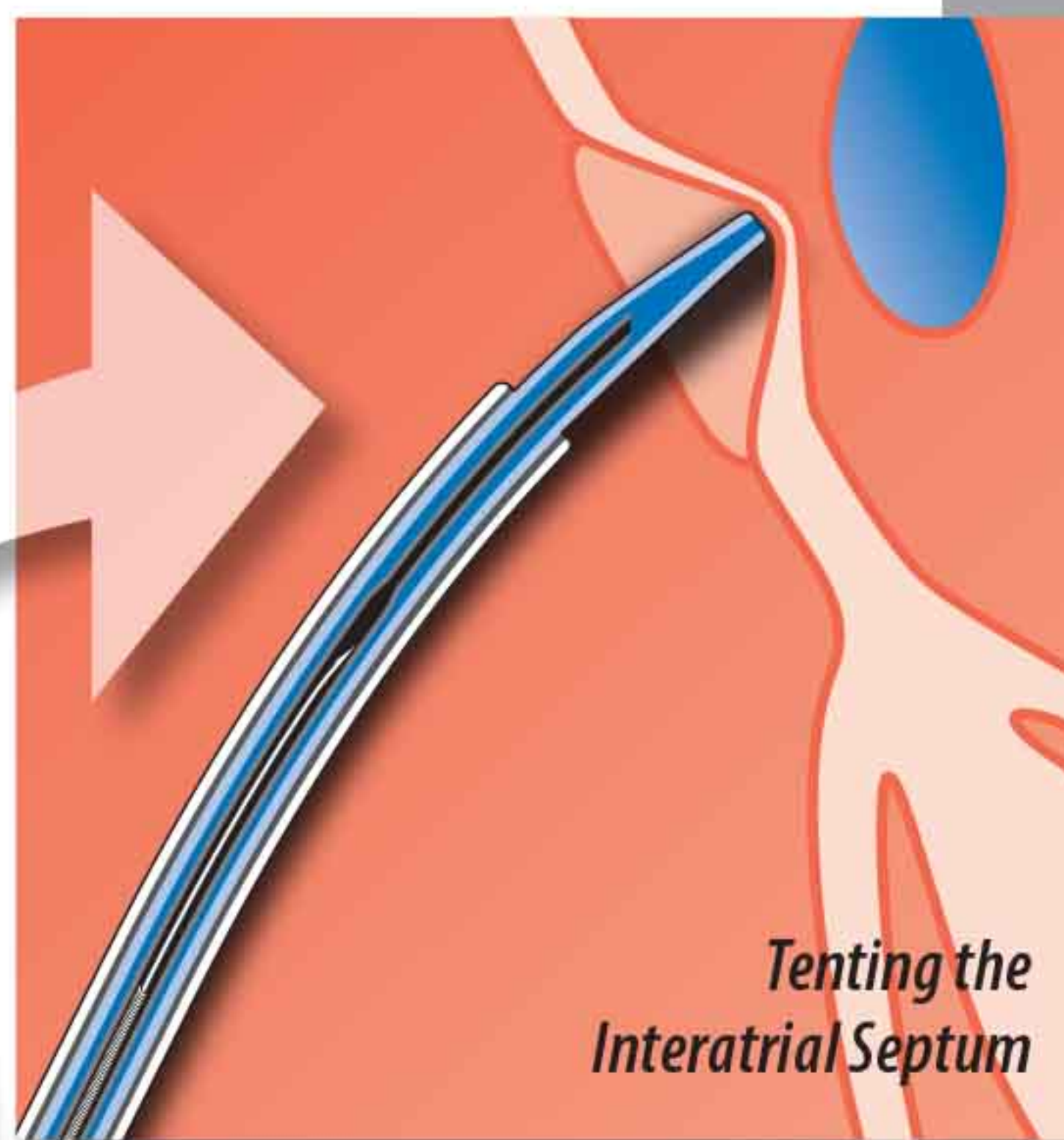

Model No.	Description	Diameter	Length
SS-140	Transseptal Guidewire with radiopaque coil	0.018"	140 cm

USA and Worldwide Patents pending. Protected by Patent Insurance.